


Faller Stylet

for

creating subcutaneous tunnel during PD Catheter insertion

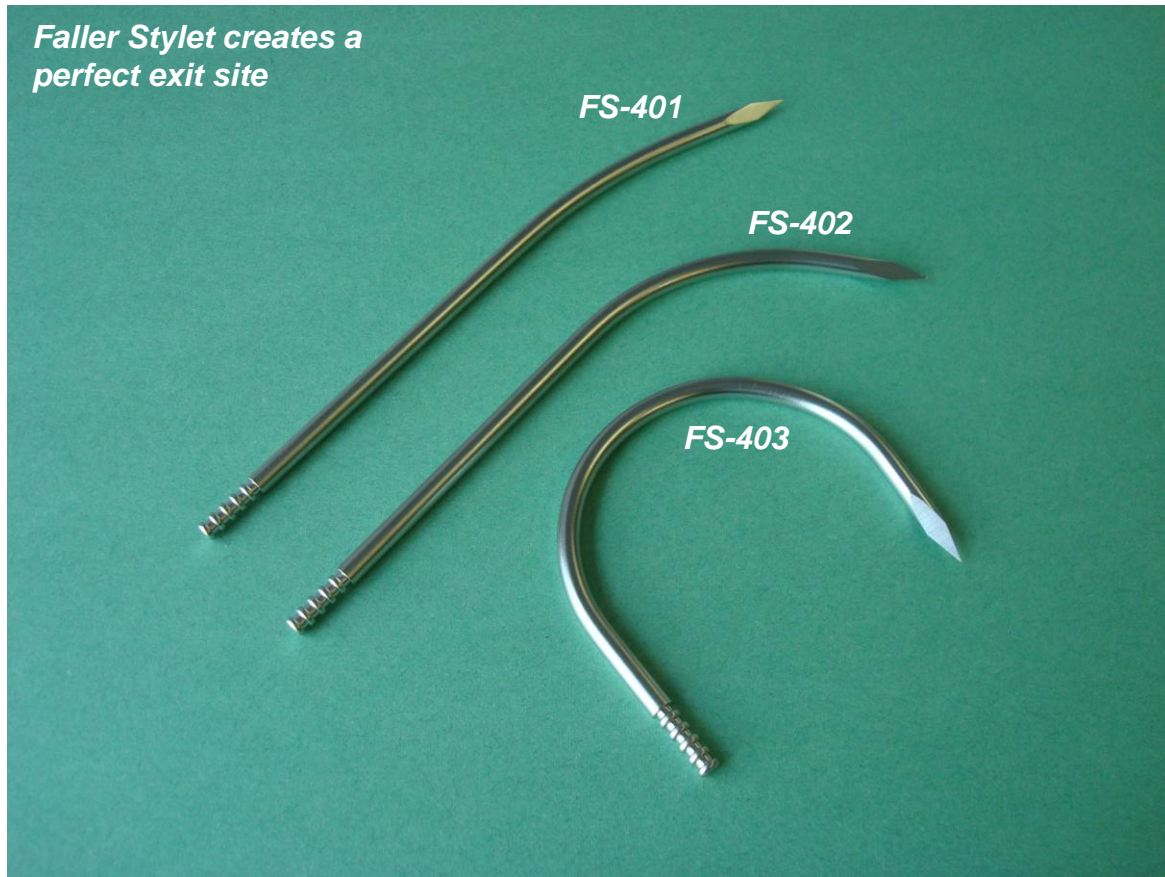

Specifications:

- *Faller stylet is a stainless steel tunneling device with a cutting edge at one end and ridges at the other end to hold the catheter;*
- *The stylet pulls the catheter through the tunnel and creates the skin exit site in one pass, in exactly the same way as a needle pulls a suture through skin;*
- *Since the faller stylet has the same outer diameter as the Medionics' PD catheter, it creates a snug exit site that requires no suturing, which minimizes exit site complications;*
- *Comes Sterile, autoclavable for possible reuse.*



Medionics International Inc., 114 Anderson Avenue, Markham, Ontario L6E 1A5 Canada
Tel: (905) 472-6544 Fax: (905) 472-6549 info@medionics.com www.medionics.com

